
## XRpad 4343 F

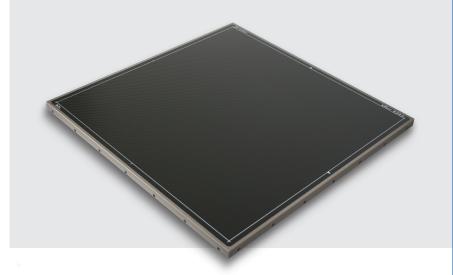

### Flat Panel X-ray Detectors

#### **Features and Benefits**

- Fixed configuration cassette detector
- True 43 cm  $\times$  43 cm (17"  $\times$  17") image
- High resolution 100 µm pixel pitch
- Up to 65,536 gray levels
- Automatic Exposure Detection (AED)
- Direct deposition CsI scintillator (XRpad 4343 F)
- High Efficiency Gadox scintillator (XRpad 4343 F-N)
- High Resolution Gadox scintillator (XRpad 4343 F-J)

#### Applications<sup>1</sup>

• Digital radiography

# Upgrade to Digital Radiography

#### **Overview**

PerkinElmer XRpad™ 4343 F is a lightweight, cassette-sized flat panel detectors for digital radiography. It fits into a conventional table or wall-stand bucky, just like a film-screen cassette.

Featuring a 18.7 million pixel image matrix and best-in-class 100  $\mu$ m pixel pitch, this detector provides exceptional image quality. Both CsI and Gadox scintillators are available. The true cassette-sized construction, with 43  $\times$  43 cm<sup>2</sup> imaging area, fits in a standard bucky, enabling easier placement and cleaning. Automatic Exposure Detection simplifies integration.

We have over 25 years of experience partnering with customers to develop products in a wide range of X-ray applications. Let our digital imaging expertise work for you.




# XRpad 4343 F

Sensor				
Panel	Single substrate amorphous silicon active TFT/diode array			
Scintillator	XRpad 4343 F: Direct Deposition CsI:TI XRpad 4343 F-N: High Efficiency Gadox XRpad 4343 F-J: High Resolution Gadox			
Pixel Matrix	4318 × 4320			
Pixel Pitch	100 µm			
Electronics				
Amplifiers	Low noise ASICs with user selectable gains			
ADC	16-bit			
Image Transfer Time	600 ms			
On-board Memory	1 GB DDR3, 4/8 GB SDHC card			
Mechanical				
Size	43 cm $\times$ 43 cm (17" $\times$ 17") cassette size			
Active Area	432 mm × 432 mm			
External Dimensions	60 mm (w) $\times$ 460 mm (l) $\times$ 15 mm (h)			
Weight	4.1 kg (9.0 lbs)			
Housing	Aluminum frame with carbon-fiber entrance window			
Communications				
Wired Data I/F	GigE via power and communication tether			
X-ray I/F	Integrated X-ray trigger control Automatic Exposure Detection			

Imaging Performance				
Typical DQE	cy/mn	XRpad 4343 F (CsI)	XRpad 4343 F-N (Gadox)	
	0	75%	35%	
	1	60%	25%	
	3	40%	8%	
Typical MTF	cy/mn	XRpad 4343 F (CsI)	XRpad 4343 F-N (Gadox)	
	1	70%	56%	
	2	40%	25%	
	4	15%	5%	
Limiting Resolution	5 cy/mm			
Energy Range	20 - 150 kV			
Environmental				
Temperature	10 – 35 °C operating			
Humidity	30 - 70 % RH operating (non-condensing)			
Accessories				
Interface and Power Unit	XRpad IPU external power supply 100 - 240 V AC, GigE, and X-ray I/F			
Regulatory				
Standards	IEC 60601-1, IEC 60601-1-2, IEC 60601-1-6, EN 62311, ISO 10993-5, ISO 10993-10			

1: Unless otherwise specified, PerkinElmer Flat Panel X-ray Detectors are components intended to be integrated into products by X-ray system manufacturers. System manufacturers are responsible for qualifying and validating their products for their intended uses and meeting all applicable regulatory requirements.

Contents in this document are subject to change without notice.

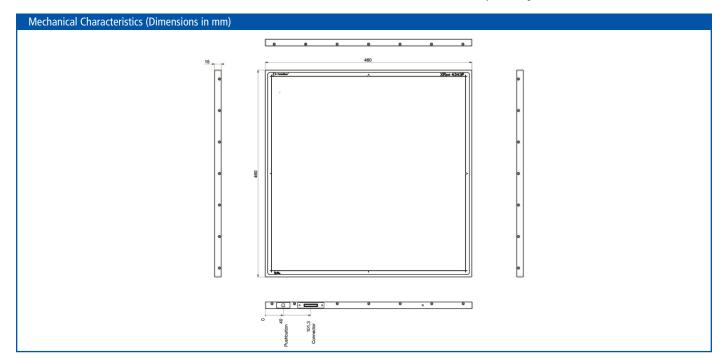

PerkinElmer Medical Imaging, LLC 2175 Mission College Blvd Santa Clara, CA 95054 USA P: +1-844-726-8228 F: +1 408-969-6493 Germany
PerkinElmer
Technologies GmbH & Co. KG
In der Rehbach 22
65396 Walluf Germany
P: +49 6123 971-300
F: +49 6123 971-600

United Kingdom Dexela Limited A PerkinElmer Company Wenlock Business Centre 50-52 Wharf Road London N1 7EU United Kingdom P: +44 20 7148 3107

China
PerkinElmer, Inc.
No. 1670 Zhangheng Road
Zhangjiang Hi-Tech Park
Shanghai 201203, PRC
P: +86 21 60645611
F: +86 21 60645666

PerkinElmer

 $Visit\ our\ website\ at\ www.perkinelmer.com/imaging-components/\ and\ contact\ us\ through\ the\ contact\ form\ provided.$